

Simulation-based prostate enucleation training: Initial experience using 3D-printed organ phantoms

Claudia Deyirmendjian¹, David-Dan Nguyen², Sero Andonian³, Melanie Aubé-Peterkin³, Julien Letendre⁴, Dean Elterman⁵, Kevin C. Zorn⁶, Bilal Chughtai⁷, Arkadiusz Miernik⁸, Andreas J. Gross⁹, Naeem Bhojani⁶

¹Faculty of Medicine, Université de Montréal, Montreal, QC, Canada; ²Division of Urology, University of Toronto, Toronto, ON, Canada; ³Division of Urology, McGill University Health Centre, Montreal, QC, Canada; ⁴Division of Urology, Maisonneuve-Rosemont Hospital, Montreal, QC, Canada; ⁵Division of Urology, University Health Network, Toronto, ON, Canada; ⁶Division of Urology, Centre Hospitalier de l'Université de Montréal, Montreal, QC, Canada; ⁷Department of Urology, Weill Cornell Medical College, New York, NY, United States; ⁸Department of Urology, University of Freiburg, Freiburg, Germany; ⁹Department of Urology, Asklepios Hospital Barmbek, Hamburg, Germany

Cite as: Deyirmendjian C, Nguyen D-D, Andonian S, et al. Simulation-based prostate enucleation training: Initial experience using 3D-printed organ phantoms. *Can Urol Assoc J* 2022;16(12):409-16. <http://dx.doi.org/10.5489/cuaj.7838>

Published online July 21, 2022

See related commentary on page 417

Abstract

Introduction: Anatomical endoscopic enucleation of the prostate (AEEP) is an effective treatment for benign prostatic hyperplasia (BPH); however, there is controversy regarding the difficulty of learning such a technique. Simulation-based training can mimic real-life surgeries and help surgeons develop skills they can transfer to the operating room, thereby improving patient safety. This study aimed to evaluate the validity of a novel organ phantom for use in AEEP simulation training.

Methods: Participants performed AEEP on organ phantom simulators during a Masterclass using one of three energy modalities: holmium:YAG laser, thulium fiber laser, or bipolar energy. The organ phantom is composed of hydrogels and uses 3D molds to recreate prostatic tissue. Participants completed a questionnaire assessing content validity, face validity, feasibility, and acceptability of using the prostate organ phantom.

Results: The novice group consisted of 13 urologists. The median number of AEEP previously performed was 0 (interquartile range [IQR] 0–2). Two experts in AEEP (surgeons having performed over 100 AEEP interventions) also participated. All participants agreed or strongly agreed that there is a role for simulators in AEEP training. Participants positively rated the overall operative experience (7.3/10). Morcellation (4.7/10) and hemostasis (3.1/10) were deemed the least realistic steps. All participants considered it feasible to incorporate this organ phantom into training programs and 92.9% agreed that it teaches skills transferrable to the operating room.

Conclusions: This study has established content and face validity for AEEP with three different energy sources for an organ phantom. Participants considered its use both feasible and appropriate for AEEP training purposes.

KEY MESSAGES

- Simulation-based training with a prostate organ phantom is feasible for novice urologists to learn AEEP.
- Content validity was established and respondents agreed there is a role for a validated simulator in AEEP training.
- Face validity was acceptable for most components of AEEP surgery (but less appropriate for hemostasis or morcellation).
- The most challenging part of enucleation using the simulator was the en-bloc technique; the second was finding the plane of enucleation.
- The 3D-printed prostate organ phantom is an accessible and reproducible simulator that 92.9% of participants believe teaches transferrable skills useful in the OR.

Introduction

Anatomical endoscopic enucleation of the prostate (AEEP) became part of the therapeutic arsenal for benign prostatic hyperplasia (BPH) about 20 years ago.^{1,2} AEEP is associated with reduced blood loss and shorter hospital stay than transurethral resection of the prostate (TURP) and open simple prostatectomy (OSP), the current respective gold standards for prostate volumes <80 cc and >80 cc, respectively.^{3,4} Although AEEP has been shown to have excellent perioperative and postoperative outcomes with outstanding long-term durability, its acceptance within the urological community has been slow.⁵⁻⁷ Some speculate that the reluctance to adopt AEEP may be due to a steep learning curve, despite some studies suggesting a similar learning curve to TURP.⁸ There is evidence that a surgeon must have completed between

30 and 50 procedures to carry out holmium AEEP safely and efficiently.⁹⁻¹¹ Other considerations affecting individual learning curve include surgeon BPH procedure volumes, endoscopic surgical skill set, and mentorship or fellowship-dedicated training.¹²⁻¹⁴

Simulation-based training has been widely proposed as a method for practicing surgeons to learn AEEP outside of fellowship tutelage.¹⁵⁻¹⁷ Integrating simulators in surgical training allows urologists to safely develop the fine-motor skills required for AEEP without adverse consequences on real patients. Different simulators have been validated for AEEP. These simulators face similar challenges of creating a realistic experience mimicking real-life AEEP and helping surgeons develop skills they can transfer to the operating room.

Among existing AEEP simulators, a prostate bench model by Matsuda and colleagues has acceptable laser-tissue interaction, irrigation and bubbles, yet is limited by lack of bleeding or morcellation.^{15,18} Virtual-reality simulators can create reproducible surgical experiences without the need for external assessment, but trainees have to search elsewhere to gain experience using real holmium lasers, irrigation fluid, or experience practicing effective laser safety protocol.¹⁹ Human cadavers are another tool for surgical training in AEEP that can create a highly realistic operative experience, including equipment, setup, and morcellation; however, they lack bleeding and their availability is limited.^{20,21}

The novel prostate organ phantom engineered by the Max Planck Institute is composed of hydrogels and uses 3D moulds to recreate prostatic tissue and anatomy. While this model has been used successfully to teach TURP to urologists, it has not yet been validated for AEEP surgical training.²² As such, this study sought to validate the 3D prostate organ phantom for AEEP training among urologists with various levels of surgical experience. In the setting of a Masterclass, participants assessed the content validity, face validity, feasibility, and acceptability of incorporating this organ phantom into AEEP training.

Methods

Study setting

An AEEP Masterclass was hosted by the Canadian Urological Association on November 23–24, 2021. Canadian urologists interested in learning AEEP were invited to participate in this accredited Masterclass. They participated in a series of lectures regarding AEEP, including equipment, pre- and postoperative care, outcomes, and advanced technique. They then practiced directly on the organ phantom simulators under the supervision of AEEP experts that have performed at least 100 AEEP cases. Each trainee was given a 30-minute timeframe to operate and was mentored by one of

the experts. The technique used was the en-bloc enucleation as described by Scoffone et al.²³ A video of the Masterclass with examples of technique can be found at <https://player.vimeo.com/video/653401942>.

Prostate organ phantom

To simulate the AEEP procedure, the “Endo Urology Trainer” was used, which is an organ phantom of the full urinary tract provided by the Max Planck Institute for Intelligent Systems and the University of Stuttgart. The 3D-printed organ phantoms are patented (EP 3251811 and WO 2017/207361) and recreate the internal anatomy of a bilobar prostate, including the verumontanum, and use two different materials to clearly distinguish the central from peripheral zones of the prostate. Images of the organ phantom before and after complete enucleation are depicted in Figure 1. The prostate phantom is designed to simulate endourological resections with the 2 µm lasers and bipolar electrocautery instruments.²²

For the purposes of this study, surgeons were assigned one of three energy modalities — holmium:YAG laser, thulium fiber laser, or bipolar energy — and performed the enucleation with genuine instruments and irrigation on the organ phantoms. For AEEP performed using the holmium and thulium fiber lasers, the Shark® continuous irrigation resectoscope with 30° telescope, 26 Fr outer sheath, and 24 Fr inner sheath with 600 µm straight distal end and 24 Fr obturator (Richard Wolf GmbH, Knittlingen, Germany) were used. The Shark® resectoscope for bipolar enucleation had a 30° telescope, 26 Fr outer sheath, 24 Fr obturator and viewing obturator with a 24 Fr bipolar enucleation electrode. The Piranha® morcellation system (Richard Wolf GmbH,

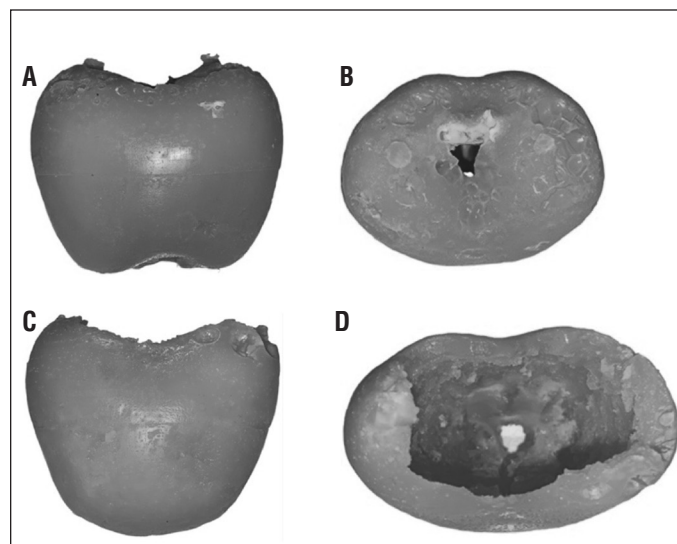


Figure 1. Organ phantom before resection (A) anterior view; (B) base view; and after complete enucleation with a holmium laser; (C) anterior view; (D) base view.

Knittlingen, Germany) was used for all stations. Figure 2 shows the setup of the simulator.

Questionnaire

Consent was obtained from each participant according to ethics board authorization (study number: 2022-10107). At the end of the second day of the Masterclass, participants completed a quantitative questionnaire evaluating the validity of the organ phantoms for use in AEEP training. The questionnaire used was a modified version of a validated questionnaire by Antunes et al.¹⁵ Demographic variables and prior experience of participants were also collected.

Content validity was measured by rating the level of agreement with the following four statements: 1) There is a role for a validated AEEP simulator in training; 2) There is a role for simulators in surgical training in general; 3) Simulation-based training and assessment for AEEP is essential to patient safety; and 4) AEEP is an effective method of treatment. Face validity was assessed by both the experts and novice group rating the realism of each component of the operative experience and each step of enucleation on a 10-point Likert scale (1 being “poorly reproduced” and 10 being “realistically reproduced”).

Feasibility was assessed posing two questions: 1) Is it feasible to incorporate this organ phantom simulator in a training program? and 2) Is it feasible to adopt the prostatic organ phantom to train and assess urologists in: a) anatomy identification; b) power settings; c) fiber positioning; d) effective technique; e) preventing injury; f) avoiding instrument damage; and g) avoiding blood loss. We assessed acceptability with the following questions: 1) Should simulation be integrated into training programs?; 2) Should simulation be part of certification and recertification?; and 3) Does this 3D model system teach transferrable skills applicable in the operating room?



Figure 2. Set-up of the simulator during the Masterclass.

Participants were asked to rate the difficulty of each step of AEEP using a five-point Likert scale.

Operative outcomes

Specimens were examined postoperatively to determine any differences between results applying the three different energy modalities. Rates of capsule perforation, as well as each specimen's approximate enucleated percentage were documented.

Statistical analysis

Statistical analysis was performed using SPSS version 27 software (SPSS, IBM Corp, Armonk, NY, U.S.). Results were summarized descriptively. Previous experience of the novice group was compared to the expert group using the Fisher's exact test. The Mann-Whitney U test was used to compare median values. The difference between mean difficulty scores was calculated with the Student's t-test. The Kruskal-Wallis test was used to determine the effect of prior experience on simulation difficulty. Statistical significance was set at a two-sided $p < 0.05$.

Results

Demographics

Thirteen participants with little or no AEEP experience (novice group) and two AEEP experts were recruited. The average age of participants was 43 years, and 100% of participants were male. Most participants were urologists with an average 11.7 years of practice, and only one participant was a fifth-year resident in urology. The median number of AEEP performed by the novice group was 0 (interquartile range [IQR] 0–2) but 12 out of 13 had performed at least 20 TURP annually, and four had performed at least 20 transurethral vaporizations annually. The vast majority (92.3%) of AEEP-novice surgeons and all experts had already used a simulator for training purposes in surgery, but this was the first experience for all participants using an AEEP-specific simulator. See Table 1 for additional details on participants' prior surgical experience.

Operative outcomes

Twenty-two prostates were enucleated: six with a holmium laser, 14 with a thulium fiber laser, and two with bipolar energy (Table 2). Organ phantoms enucleated with a holmium laser revealed the fewest capsule perforations, with only 16.7% of specimens being perforated, compared to 57.1% of thulium fiber laser enucleations and all bipolar enucleations.

Table 1. Previous experience of participants in the Masterclass

Previous experience		Expertise		p
		Novice (n=13)	Experts (n=2)	
Duration of urological practice, years (SD)		12.9 (11.2)	5.0 (4.2)	0.358
Number of TURP procedures performed annually (%)	Less than 20	1 (7.7)	2 (100)	0.038
	20–50	10 (76.9)	0	
	More than 50	2 (15.4)	0	
Number of transurethral vaporisation procedures performed annually (%)	Less than 20	9 (69.2)	2 (100)	1.000
	20–50	4 (30.8)	0	
	More than 50	0	0	
Percentage of practice that consisted of laser enucleation of the prostate procedures in a year (%)	Less than 50%	13 (100)	0	0.010
	About 50%	0	1 (50)	
	More than 50%	0	1 (50)	
Median approximate number of laser enucleation of the prostate procedures performed (IQR)		0 (0–2)	560 (NA)	0.010
Median number of years of experience practicing laser enucleation of the prostate (IQR)		0 (0–0.5)	5 (NA)	0.019
Type of affiliated institution	Non-teaching	6 (46.2)	0	0.486
	Teaching	7 (53.8)	2 (100)	
Before this Masterclass, have you ever used a simulator for surgical training?	No	1 (7.7)	0	1.000
	Yes	12 (92.3)	2 (100)	
Before this Masterclass, have you ever used a simulator for laser enucleation of the prostate training?	No	13 (100)	2 (100)	1.000
	Yes	0	0	

Boldface print indicates significance. IQR: interquartile range; SD: standard deviation; TURP: transurethral resection of the prostate.

Content validity

Regarding content validity, 100% of participants agreed or strongly agreed that AEEP was an effective treatment for BPH (Figure 3). All participants also agreed or strongly agreed that simulators have a role to play in surgical training, and specifically in AEEP training. Most agreed that simulation-based training is essential to patient safety, but 13.3% disagreed with this statement.

Face validity

Face validity was acceptable for the overall operative experience, which was rated 7.3 on a 10-point scale (Figure 4A). Instrumentation was the most realistic component (9.6/10), while laser-tissue interaction obtained the lowest score

(6.1/10). The realism of each individual enucleation step was also rated; creating the 5, 6, or 7 o'clock groove and the anterior groove were judged the most realistic (7.9/10) (Figure 4B). Most steps scored at least 6/10; however, hemostasis (3.1/10) and morcellation (4.7/10) were found to be poorly reproduced.

Acceptability and feasibility

It was considered acceptable by 100% of our participants that simulation be integrated into training programs (Figure 5). Only a minority of subjects (33.3%) found it appropriate to make simulation mandatory for certification and recertification. All subjects agreed that it is feasible to incorporate this organ phantom simulator within a training program, and 86.7% of participants believed it was feasible to use the simulator to train and assess urologists in identifying anatomy and effective technique. Most disagreed that the simulator was a feasible method for training blood-loss prevention (66.7%). Over 90% of participants claimed that the simulator teaches transferable skills useful in the operating room.

Difficulty

The en-bloc technique appears to be one of the most challenging enucleation steps: about three-quarters of subjects found it at least “difficult,” and it received a mean score

Table 2. Performance results of the Masterclass

Outcomes of Masterclass		Number of specimens (n=22)
Perforated specimens (% per laser type)	Holmium (n=6)	1 (16.7)
	Thulium fiber (n=14)	8 (57.1)
	Bipolar (n=2)	2 (100)
Amount resected (% of total)	Less than 50%	9 (40.9)
	About 50%	7 (31.8)
	More than 50%	6 (27.3)

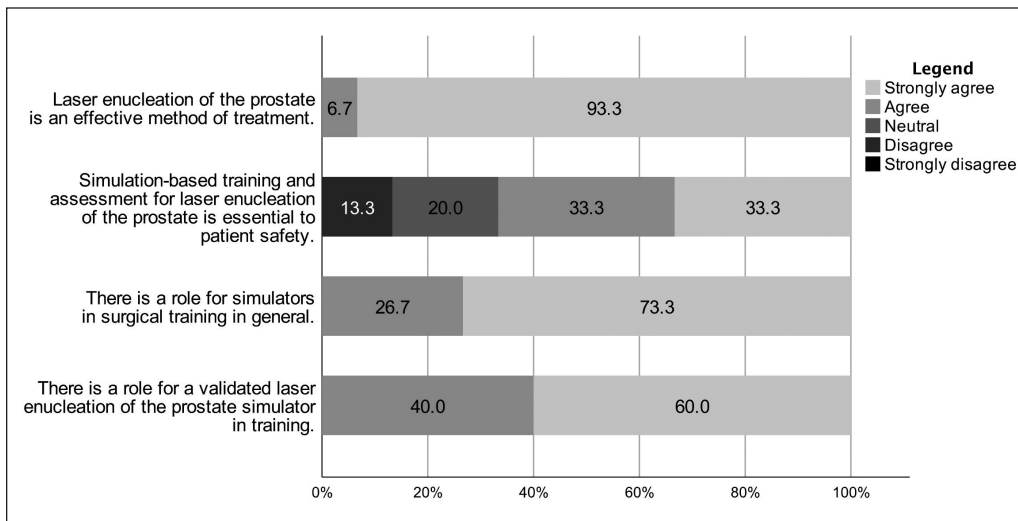


Figure 3. Content validity of the simulator.

of 3.4 on a five-point scale (range 2–5). Morcellation was also considered “extremely difficult” by 26.7% of subjects; however, many comments emphasized that the “morcellator was not working well” or “was blocked.” Three participants marked hemostasis as non-applicable. Number of TURP or vaporizations performed, years of practice, or type of teaching institution revealed no impact on difficulty. Overall difficulty scores for each step of enucleation are summarized in Table 3.

Discussion

Bioengineered 3D organ phantoms are useful tools for various applications in urology surgical education.²⁴ Notably, they can serve as patient education models, patient-specific models for rehearsing complex surgeries, and general

for renal diseases or prostate cancer; they are lacking for BPH surgical training.^{24,27}

The novel organ phantom in the present study represents important progress in training AEEP-novice urologists in an era of rapidly growing and increasingly difficult surgical techniques. In the present study, we validated the use of a novel 3D organ phantom for simulation-based AEEP training and found its use both feasible and acceptable.

Face validity involved evaluating how realistic the simulation was in terms of the operative experience and specifically in each step of enucleation. Our results are in line with findings from another biosynthetic bench model patented by Matsuda et al, which was validated by two distinct studies.²⁸ Antunes et al studied 40 urologists taking a course for holmium laser enucleation of the prostate (HoLEP), which was carried out on a trilobar prostate model.¹⁵ Content and face

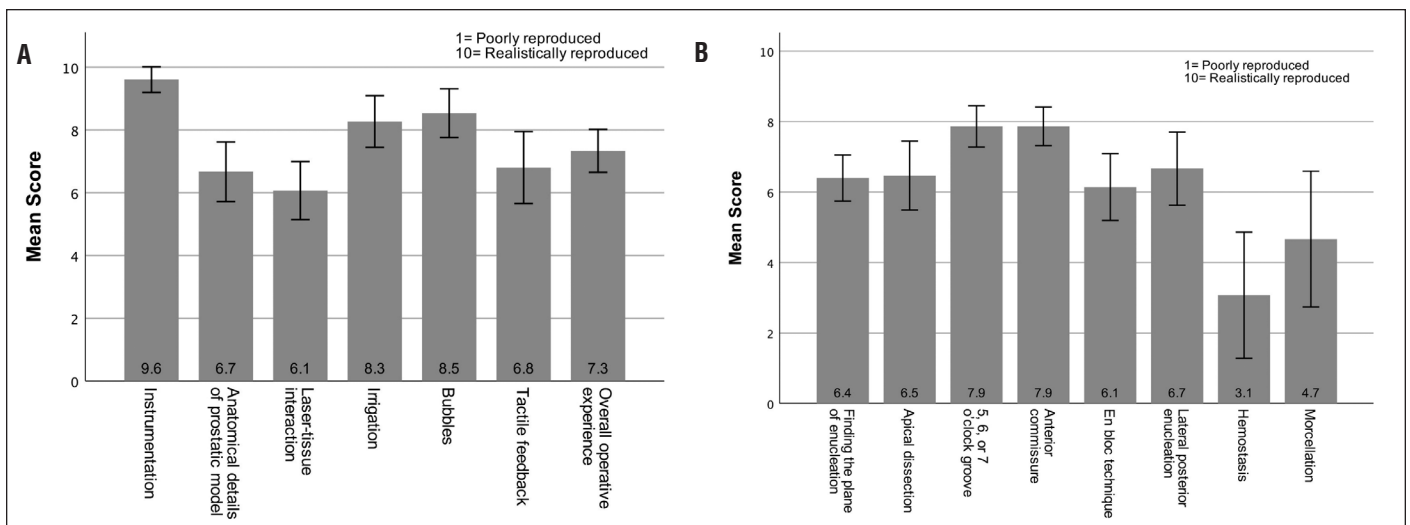


Figure 4. Face validity of (A) the operative experience and (B) each step of enucleation using the simulator.

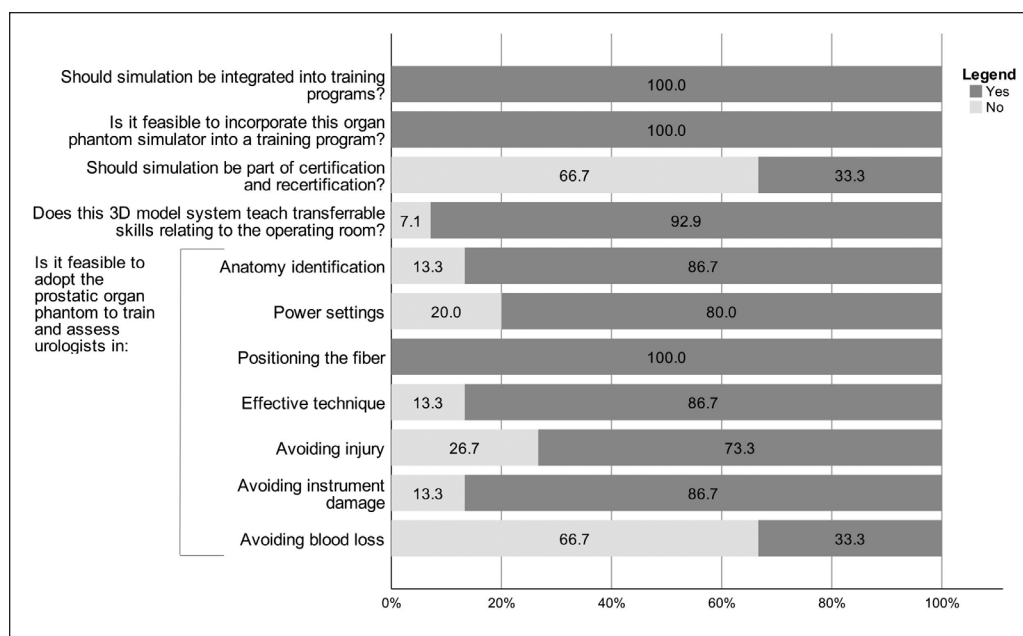


Figure 5. Feasibility and acceptability of the organ phantom in training.

validity were confirmed, but the overall operative experience was scored slightly higher (8.4 on a 10-point scale) than in our study. Instrumentation was the most realistic component, given that real equipment was used, a finding corroborated by our study. Laser-tissue interaction came in second (8.6/10), whereas laser-tissue interaction was our study's least realistic operative element (6.1/10). Note that our study featured three different energy modalities rather than just one, which may be a confounding factor.

Aydin et al conducted a prospective study of 36 urologists and residents practicing HoLEP on the same prostate model.¹⁸ Their overall operative experience was rated 7.4 on 10, which closely resembles our study finding (7.3/10). The vast majority of their participants agreed on the feasibility of applying the simulator for anatomy identification, power settings, positioning the fiber, effective technique, avoiding injury, avoiding instrument damage, and avoiding blood loss. Those findings contrast with our results, which show that two-thirds of participants failed to agree that avoiding blood loss was feasible; however, all other components were considered feasible.

Our participants confirmed the content validity of the organ phantom, although they expressed reservations about the need for simulation-based training for patient safety. Few can argue about the general utility of simulators in surgical training, but their precise role remains uncertain. This may be attributed to the fact that simulators are considered to be one adjunct among a variety of AEEP training methods.²⁹ When urologists were surveyed in the Aydin et al study about what they thought was the ideal training method, only 13% stated simulation alone and 87% believed it was supervised simulation together with operative room training.¹⁸ Other

impactful resources for learning AEEP are mentorship and proctorship, which are known to significantly lower the learning curve and reveal a positive safety profile.^{9,12,30} One study by Netsch et al found that a mentor-based approach could help urologists adopt thulium vapoenucleation and overcome possible complications at the beginning of the learning curve.¹⁴ Structured training programs can use both mentorship and simulation training, along with strategic selection of beginner cases, to lower the learning curve. One example of a structured training program for AEEP is the Holmium User Group developed in the U.K.³¹

Their approach uses modular progression to learn each step of HoLEP. Future studies should focus on longitudinal outcomes in practice when integrating both a simulator and a mentorship-based model.

Participants found the en-bloc technique to be the most challenging part of AEEP. In a training program teaching the en-bloc technique, it would thus be important to allocate enough time to learn this technique or to begin with teaching a two-lobe or three-lobe technique, as the least difficult steps are the anterior commissure and bladder neck incisions at 5, 6 or 7 o'clock.^{15,31}

Three different energy sources were accessed in this study. While AEEP-novice urologists carried out the procedures — making a higher rate of capsule perforation likely — specimens enucleated with the thulium fiber laser

Table 3. Difficulty of each step of enucleation using the simulator

Steps of enucleation	Mean score* (range)	Novice	Experts	p
Finding the plane of enucleation	2.9 (2–4)	2.8 (2–4)	3.5 (3–4)	0.207
Apical dissection	2.6 (1–4)	2.6 (1–4)	2.5 (2–3)	0.903
5, 6, or 7 o'clock groove	1.7 (1–3)	1.7 (1–3)	1.5 (1–2)	0.740
Anterior commissure	1.7 (1–4)	1.7 (1–4)	2.0 (1–3)	0.689
En bloc technique	3.4 (2–5)	3.5 (2–5)	3.0 (2–4)	0.569
Lateral posterior enucleation	2.3 (1–4)	2.2 (1–3)	2.5 (1–4)	0.887
Hemostasis	1.8 (1–5)	1.9 (1–5)	1.0 (1–1)	0.418
Morcellation	2.6 (1–5)	2.9 (1–5)	1.0 (1–1)	0.002

*On a 5-point scale, where: 1=not difficult, 2=slightly difficult, 3=difficult, 4=very difficult, 5=extremely difficult. Boldface print indicates significance.

and bipolar energy were disproportionately perforated. On the other hand, only one of six specimens enucleated by novice urologists with the holmium laser showed a capsule perforation. The thulium fiber laser is appreciated for its continuous laser, which enables easier plane correction and excellent hemostasis, and has a learning curve that may be shorter than or equivalent to HoLEP, as described by Enikeev et al.¹¹ In their randomized trial assessing the learning curve of three different energy sources with 30 patients in each group, the complication rates and specifically capsule perforation rates were very similar between the holmium:YAG and thulium fiber laser groups. Our findings may reflect the organ phantom's possible compatibility with the holmium:YAG laser since the model was developed for use with the holmium:YAG laser.

One strength of the organ phantom model in the current study was its acceptability and feasibility. The model was created via 3D printing and is easily reproducible. All participants considered it feasible to incorporate the organ phantom model into AEEP training. The organ phantom was appreciated for being particularly useful in assessments regarding: positioning the fiber, identifying anatomy, and avoiding instrument damage. The vast majority (92.9%) acknowledged that it teaches transferrable skills needed in the operating room.

Limitations

The present study is not without limitations. First, the small sample size must be considered. The Masterclass had a total of 15 participants, only two of whom were AEEP experts. Because of the few AEEP experts present, accuracy of face validity may be limited, and we were unable to assess construct validity, which should be evaluated in future research with an appropriate sample size.

Second, this prostate model was designed for HoLEP, but ultimately there were three energy modalities used in the Masterclass to allow AEEP-novice urologists to practice on different modalities. While it may have been interesting for participants to try different energy sources during the Masterclass, it is difficult to detect trends within each group due to the small sample sizes and the limited number of sessions practicing with each modality. Future studies may benefit from limiting their scope to one energy source and allowing participants more sessions to practice.

Third, while the 3D prostate organ phantom realistically reproduced prostatic anatomy, it lacks hemodynamic factors. The surgeon must rely on another resource to learn coagulation techniques.

Finally, the morcellator was not entirely functional and consequently not highly ranked for face validity. Similar criticisms regarding hemostasis and morcellation have been noted in conjunction with other synthetic models.¹⁵

Nevertheless, the present study offers insight on the first-time use of a 3D-printed organ phantom in AEEP training. This information can be relied upon when adopting these models in training programs or longitudinal training studies.

Conclusions

The 3D-bioprinted prostate organ phantom is an accessible and reproducible model that allows AEEP-novice urologists to practice this surgical technique safely. This study has established content and face validity for AEEP using the holmium:YAG laser, thulium fiber laser, and bipolar energy on a novel organ phantom. Participants considered its use both feasible and acceptable for AEEP training purposes.

Competing interests: The authors do not report any competing personal or financial interests related to this work.

Acknowledgments: The organ phantoms were engineered at the Max Planck Institute for Intelligent Systems in Stuttgart, Germany and at the University of Stuttgart. The two-day AEEP Masterclass was hosted by the Canadian Urological Association.

This paper has been peer-reviewed.

References

1. Hiraoka Y. A new method of prostatectomy, transurethral detachment, and resection of benign prostatic hyperplasia. *Nihon Ika Daigaku Zasshi* 1983;50:896-8. <https://doi.org/10.1272/jnms1923.50.896>
2. Gilling PJ, Kennett K, Das AK, et al. Holmium laser enucleation of the prostate (HoLEP) combined with transurethral tissue morcellation: An update on the early clinical experience. *J Endourol* 1998;12:457-9. <https://doi.org/10.1089/end.1998.12.457>
3. Tan A, Liao C, Mo Z, et al. Meta-analysis of holmium laser enucleation vs. transurethral resection of the prostate for symptomatic prostatic obstruction. *Br J Surg* 2007;94:1201-8. <https://doi.org/10.1002/bjs.5916>
4. Zhong J, Feng Z, Peng Y, et al. A systematic review and meta-analysis of efficacy and safety following holmium laser enucleation of prostate and transurethral resection of prostate for benign prostatic hyperplasia. *Urology* 2019;131:14-20. <https://doi.org/10.1016/j.urol.2019.03.034>
5. Robles J, Pais V, Miller N. Mind the gaps: Adoption and underutilization of holmium laser enucleation of the prostate in the United States from 2008–2014. *J Endourol* 2020;34:770-6. <https://doi.org/10.1089/end.2019.0603>
6. Wroclawski ML, Teles SB, Carneiro A. Anatomical endoscopic enucleation of the prostate: The next gold standard? No! (or not yet!). *Andrologia* 2020;52:e13707. <https://doi.org/10.1111/and.13707>
7. Ibrahim A, Alharbi M, Elhilali MM, et al. 18 years of holmium laser enucleation of the prostate: A single-center experience. *J Urol* 2019;202:795-800. <https://doi.org/10.1097/JU.000000000000280>
8. Robert G, Cornu JN, Fourmarié M, et al. Multicentre prospective evaluation of the learning curve of holmium laser enucleation of the prostate (HoLEP). *BJU Int* 2016;117:495-9. <https://doi.org/10.1111/bju.13124>
9. Kampantais S, Dimopoulos P, Tasleem A, et al. Assessing the learning curve of holmium laser enucleation of prostate (HoLEP). A systematic review. *Urology* 2018;120:9-22. <https://doi.org/10.1016/j.urol.2018.06.012>
10. Hirasawa Y, Kato Y, Fujita K. Transurethral enucleation with bipolar for benign prostatic hyperplasia: 2-year outcomes and the learning curve of a single surgeon's experience of 603 consecutive patients. *J Endourol* 2017;31:679-85. <https://doi.org/10.1089/end.2017.0092>

11. Enikeev D, Glybochko P, Rapoport L, et al. A randomized trial comparing the learning curve of 3 endoscopic enucleation techniques (HoLEP, ThuLEP, and MEP) for BPH using mentoring approach – initial results. *Urology* 2018;121:51-7. <https://doi.org/10.1016/j.urology.2018.06.045>
12. Elzayat EA, Elhilali MM. Holmium laser enucleation of the prostate (HoLEP): Long-term results, reoperation rate, and possible impact of the learning curve. *Eur Urol* 2007;52:1465-72. <https://doi.org/10.1016/j.eururo.2007.04.074>
13. Placer J, Gelabert-Mas A, Vallmanya F, et al. Holmium laser enucleation of prostate: Outcome and complications of self-taught learning curve. *Urology* 2009;73:1042-8. <https://doi.org/10.1016/j.urology.2008.12.052>
14. Netsch C, Bach T, Herrmann TRW, et al. Evaluation of the learning curve for thulium vapoenucleation of the prostate (ThuVEP) using a mentor-based approach. *World J Urol* 2013;31:1231-8. <https://doi.org/10.1007/s00345-012-0894-1>
15. Antunes AA, Iscaife A, Barbosa JABA, et al. Holmium laser enucleation of the prostate simulation: Analysis of realism and level of difficulty by holmium laser enucleation of the prostate-naïve urologists. *Urology* 2019;125:34-9. <https://doi.org/10.1016/j.urology.2018.10.055>
16. Wright HC, Fedrigon D, De S. Learning From Those who learned: A survey of fellowship-trained HoLEP surgeons and their current practice patterns. *Urology* 2021;149:193-8. <https://doi.org/10.1016/j.urology.2020.12.023>
17. Brunkhorst O, Ahmed K, Nehikare O, et al. Evaluation of the learning curve for holmium laser enucleation of the prostate using multiple outcome measures. *Urology* 2015;86:824-9. <https://doi.org/10.1016/j.urology.2015.07.021>
18. Aydin A, Ahmed K, Brewin J, et al. Face and content validation of the prostatic hyperplasia model and holmium laser surgery simulator. *J Surg Educ* 2014;71:339-44. <https://doi.org/10.1016/j.jsurg.2013.11.004>
19. Kuronen-Stewart C, Ahmed K, Aydin A, et al. Holmium laser enucleation of the prostate: Simulation-based training curriculum and validation. *Urology* 2015;86:639-46. <https://doi.org/10.1016/j.urology.2015.06.008>
20. Tunc L, Bozzini G, Scoffone CM, et al. Determination of face and content validity of cadaveric model for holmium anatomic endoscopic enucleation of the prostate training: An ESUT AEEP group study. *Eur Urol Open Sci* 2021;32:28-34. <https://doi.org/10.1016/j.euro.2021.07.005>
21. Fassi-Fehri H, Matillon X, Ye H, et al. Endoscopic holmium laser enucleation of the prostate (HoLEP) in the human cadaver: Realism and feasibility. *Clin Surg*. Published online 2021. Available at : <http://clinicsinsurgery.com/>. Accessed February 2022
22. Choi E, Adams F, Palagi S, et al. A high-fidelity phantom for the simulation and quantitative evaluation of transurethral resection of the prostate. *Ann Biomed Eng* 2020;48:437-46. <https://doi.org/10.1007/s10439-019-02361-7>
23. Scoffone CM, Cracco CM. The en-bloc no-touch holmium laser enucleation of the prostate (HoLEP) technique. *World J Urol* 2016;34:1175-81. <https://doi.org/10.1007/s00345-015-1741-y>
24. Agung NP, Nadhif MH, Irdam GA, et al. The role of 3D-printed phantoms and devices for organ-specified appliances in urology. *Int J Bioprinting* 2021;7:1-14. <https://doi.org/10.18063/ijb.v7i2.333>
25. Ghazi AE, Teplitz BA. Role of 3D printing in surgical education for robotic urology procedures. *Transl Androl Urol* 2020;9:931-41. <https://doi.org/10.21037/tau.2020.01.03>
26. Chen MY, Skewes J, Desselle M, et al. Current applications of three-dimensional printing in urology. *BJU Int* 2020;125:17-27. <https://doi.org/10.1111/bju.14928>
27. Tatar İ, Huri E, Selçuk İ, et al. Review of the effect of 3d medical printing and virtual reality on urology training with 'MedRain3DModsim' Erasmus + European Union Project. *Turkish J Med Sci* 2019;49:1257-70. <https://doi.org/10.3906/sag-1905-73>
28. Matsuda K, Kinoshita H, Okamoto Y. Prostatic hyperplasia model and prostate surgery simulator. Published online May 23, 2013. Available at : <https://patentscope.wipo.int/search/en/detail.jsf?docId=W02013073210&tab=PCTBIBLIO>. Accessed July 8, 2021
29. Teoh JYC, Cho CL, Wei Y, et al. Surgical training for anatomical endoscopic enucleation of the prostate. *Andrologia* 2020;52:e13708. <https://doi.org/10.1111/and.13708>
30. Westhofen T, Weinhold P, Kolb M, et al. Evaluation of holmium laser enucleation of the prostate learning curves with and without a structured training programme. *Curr Urol* 2020;14:191-9. <https://doi.org/10.1159/000499239>
31. Aho T, Herrmann TRW. Description of a modular mentorship program for holmium laser enucleation of the prostate. *World J Urol* 2015;33:497-502. <https://doi.org/10.1007/s00345-014-1407-1>

Correspondence: Dr. Naeem Bhojani, Division of Urology, Centre Hospitalier de l'Université de Montréal, Montreal, QC, Canada; naeem.bhojani@gmail.com